

clinell®

www.clinell.com.au

SUPER
ABSORBENT
SOAKS UP TO
1 LITRE OF
LIQUID



SOAKS UP SPILLS SAFELY, IN SECONDS

Proven to kill microorganisms associated with blood, urine and faecal spills

Including Hepatitis B & C, HIV, *Clostridium difficile*, Norovirus and many more. High power oxidative kill against all microorganisms including non-enveloped viruses and bacterial endospores.

Simple, quick and highly effective

Replacing traditional chlorine based spill kits which are time consuming and hazardous to use. Due to its exceptional ease of use it increases compliance which reduces incidents of pathogen transmission.

Easy and safe disposal of soiled wipes

Clinell Spill Wipes packaging is designed so that all items, including the protective clothing worn by the user, can be disposed of in the re-sealable bag and placed into the appropriate disposal channel.

Based on technology from our leading Sporicidal Wipes

Similar to Clinell Sporicidal Wipes, the super absorbent pad generates peracetic acid and hydrogen peroxide on activation. This proven formula is also safe to use on urine, unlike traditional chlorine based products.

SPILL WIPES

Clinell Spill Wipes are designed to make dealing with body fluid spills quick and easy. The pack contains one super absorbent, peracetic acid generating pad and two large disinfectant wipes. All within a re-sealable bag for the easy and safe disposal of waste.



The traditional body fluid spills protocol requires the spill to be decontaminated with chlorine before the spill can be safely removed. This involves a lengthy procedure; place a barrier agent around the edges of the spill, pour 10,000ppm of chlorine granules over the spill and leave for a few minutes to neutralise the

blood-borne pathogens. Next, remove the chlorine granules with a scooper/ dustpan and dispose of in a biohazard bag. Wipe the area with additional absorbent towels to remove any residue.

This entire process takes at least 10 minutes and exposes the user to

hazardous chemicals as well as carcinogenic toxic fumes. This may be why compliance is very low with current chlorine blood spill kits. In addition, chlorine spill kits cannot be used on urine - meaning two different spill kits for blood and urine are necessary. With Clinell Spill Wipes you only need one.

Clinell Spill Wipes were created to make the entire process of dealing with body fluid spills easier and safer. They increase compliance in this necessary, yet often neglected, infection control procedure in a cost effective manner.



1. Wear appropriate PPE.



2. Tear the pack open.




3. Remove the wipes.



4. A is the active side and should be placed on the spill.



5. Place the active side (A) face down onto the spill. Leave to absorb for 30 secs.

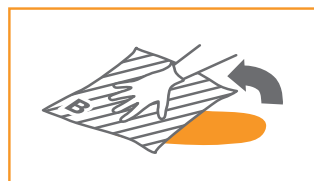

6. Push down on the plastic backed side (B) and move around until the spill is fully absorbed.



7. Place the soiled wipe back into the pack.



8. Remove a disinfectant wipe from the sachet.



9. Clean the spill area in an 'S' shaped motion, from clean to dirty.



10. Place the soiled wipe and empty sachet back into the pack.



11. If required repeat steps 8-10 with the second wipe and re-seal.



12. Dispose of in the appropriate waste stream.

PRODUCT INFORMATION

Spill Wipes

Wall Mounted Dispenser

UNIT OF ISSUE

Single unit

(for CSW1AUS)

ORDER CODE

CSW1AUS

CSW1D

GHA190019