



ProphylCare®

Directa's products are developed, tested and evaluated by our selected team of highly-qualified and renowned Swedish and international dental professionals. The aim is to make life easier for all dental professionals in their daily clinical work.


DIRECTA



The perfect starter, effective removal of stubborn stains

For discoloured surfaces and thick layers of plaque

Prod No	Product Description
731113	PropyCare blue RDA 250, 125 micron, tubes, 96 g 60 ml
690123	PropyCare Single Dose, blue 250, 144x 2 g
690119	PropyCare Single Dose, mix, 144x 2 g

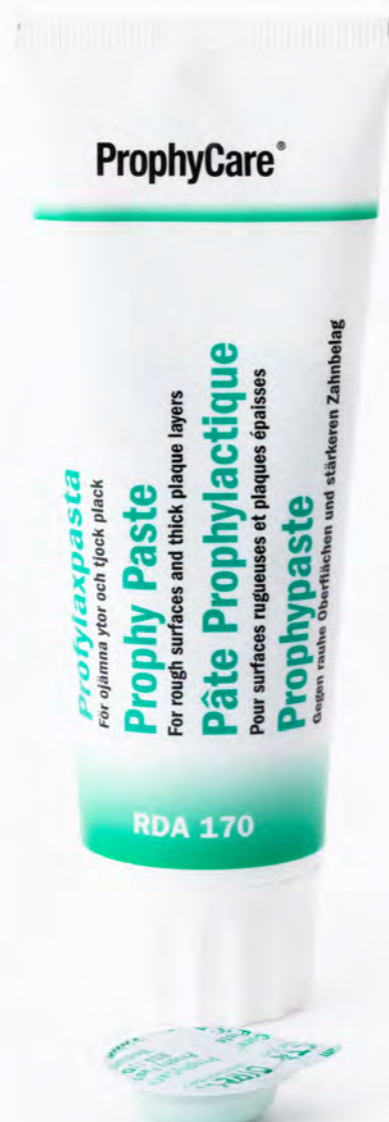

Less abrasive, but still effective

For thick layers of plaque and discoloured surfaces

Prod No	Product Description
731114	PropyCare green RDA 170, 70 micron, tubes, 97 g 60 ml
690122	PropyCare Single Dose, green 170, 144x 2 g
690119	PropyCare Single Dose, mix, 144x 2 g



Gentle cleaning, tissue friendly stain remover

For small amounts of plaque, root surfaces and lightly discoloured teeth

Prod No	Product Description
731112	PropyCare red RDA 120, 20 micron, tubes, 84 g 60 ml
690121	PropyCare Single Dose, red 120, 144x 2 g
690119	PropyCare Single Dose, mix, 144x 2 g



Ideal for final polishing

For final polishing of teeth and fillings

Prod No	Product Description
731111	PropyCare yellow RDA 40, 2 micron, tubes, 84 g 60 ml
690120	PropyCare Single Dose, yellow 40, 144x 2 g
690119	PropyCare Single Dose, mix, 144x 2 g

ProphyCare® Single Dose

ProphyCare Single Dose is available in assortment or refill packages. The box is designed to keep the single dose cups tidily in place and give an overview of the contents, allowing users to quickly select what is needed for each treatment. Choose ProphyCare after evaluating the degree of discolouration, presence of exposed dentin or cementum, aesthetic and other restorative materials, areas of demineralisation, etc. Use at least two different ProphyCare RDA values for each treatment - RDA 250 and/or RDA 170 for cleaning and RDA 120 and/or RDA 40 for polishing.

Prod No	Product Description
690123	ProphyCare Blue 250, single dose, 144x 2 g
690122	ProphyCare Green 170, single dose, 144x 2 g
690121	ProphyCare Red 120, single dose, 144x 2 g
690120	ProphyCare Yellow 40, single dose, 144x 2 g
690119	ProphyCare Mix, single dose, 144x 2 g
690125	ProphyCare PRO single Dose, 144x 2 g

Extremely hygienic

You are always using a freshly sealed package. No risk of cross-contamination, ensuring high hygiene standards for your patient.

Adhesive base

The single dose cup has an adhesive base that easily sticks to the glove or any other surface.

No waste

The shallow packaging with rounded edges is designed to minimize waste. All the content can easily be collected by the polishing cup.

Perfect size

Each package contains the ideal quantity for one treatment.




Selection Guide

	Blue RDA 250	Green RDA 170	Red RDA 120	Yellow RDA 40
Starter	●	●		
Very discoloured surfaces	●			
Thick plaque	●	●		
Discoloured teeth		●		
Small amounts of plaque			●	
Root surfaces			●	
Slightly discoloured teeth			●	
Fillings				●
Pediatric				●
Cervical toothbrush abrasion lesions				●
Implants			●	●
Final polishing			●	●
After air polishing treatment			●	●
Final polish				●

ProphyCare® Blue 250
Blue Coarse, Paraben Free
RDA 250, 125 micron

Perfect starter
For discoloured surfaces and thick layers of plaque
Effective removal of stubborn stains caused by agents such as coffee and tobacco
Abrasive agent: pumice
Contains 0,1% fluoride
Use in combination with other ProphyCare paste

Prod No	Product Description
731116	ProphyCare Blue 250, Paraben Free, tubes, 96 g 60 ml

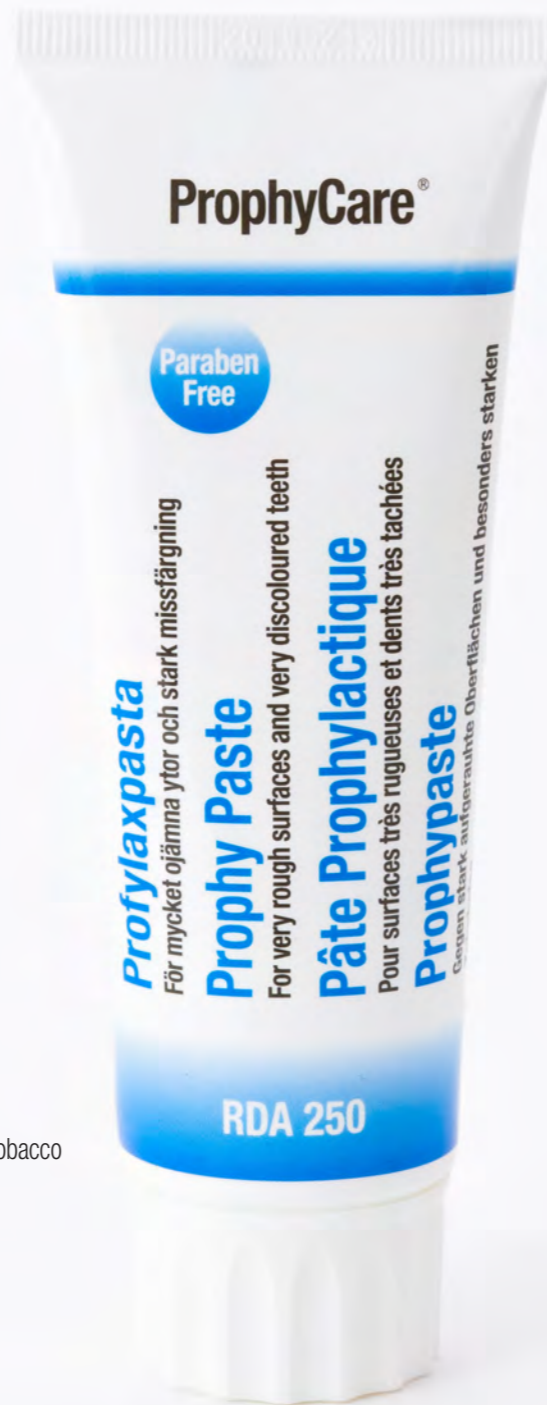

ProphyCare® Red 120
Red Fine, Paraben Free
RDA 120, 20 micron

For removing plaque and surface discolourations
For small amount of plaque, root surfaces and lightly discoloured teeth
Excellent for polishing
Abrasive agent: silica
Contains 0,1% fluoride
Use in combination with other ProphyCare Paste

Prod No	Product Description
731115	ProphyCare Red 120, Paraben Free, tubes, 84 g 60 ml



Desensitisation

ProphyCare HAp contains hydroxyapatite that effectively contributes to relief from hypersensitivity.

Remineralising

Studies* have shown that products containing hydroxyapatite show better remineralising properties than traditional fluoride products.

Smooth & Creamy

The smooth and creamy texture is easy to rinse off and does not splatter.

One Step Procedure

The same paste is used for effective stain removal and gentle final polishing – in one step. ProphyCare HAp's high cleaning power converts into gentle polishing as the sizes of the particles are fragmented.

*Enamel and dentine remineralization by nano-hydroxyapatite toothpastes, JOD, Volume 39, Issue 6, June 2011, Pages 430-437, study by Peter Tshoppe, Daniela L.Zandim, Peter Martus, Andrej M.Kielbassa

*Randomized, placebo-controlled study of the efficacy of a calcium phosphate containing paste on dentin hypersensitivity, Dent Mater. 2015 Nov;31(11):1298-303. doi: 10.1016/j.dental.2015.08.162. Epub 2015 Sep 8. Study by Deepak Mehta, Vishwas Gowda, Werner J Finger, Keiichi Sasaki

ProphyCare® HAp	
Starter	●
Very discoloured surfaces	
Thick plaque	
Discoloured teeth	●
Universal paste	●
Small amounts of plaque	●
Root surfaces	●
Slightly discoloured teeth	●
Fillings	●
Pediatric	●
Cervical toothbrush abrasion lesions	●
Implants	●
Final polishing	●
After air polishing treatment	●
To be used for polishing after scaling to relieve hypersensitivity after treatment	●

ProphyCare® HAp

With ProphyCare HAp we are offering the new generation of prophy pastes, providing your patient with a powerful relief from dentin hypersensitivity.

The vision was to create a new generation of prophy paste that combines all the features you expect from your paste, with a proven capacity to relieve dentin hypersensitivity.

ProphyCare HAp has a unique composition including hydroxyapatite that effectively will block the open dentin tubules on an exposed dentin surface.

The result is a barrier that will provide the patient with instant relief from hypersensitivity.

ProphyCare HAp is a one-step prophy paste for outstanding cleaning and polishing results.

Prod No	Product Description
731117	ProphyCare HAp Tubes, 60g 60 ml



Excellent cleaning effect with minimal abrasion

ProphyCare PRO offers effective stain removal properties at the beginning of treatment and decreases progressively to a lower RDA value, giving the paste its final polishing properties at the end of treatment.

One step procedure

ProphyCare PRO is the one-step prophy paste used for achieving outstanding cleaning and polishing results, with reduced treatment time. The same paste is used for both stain removal and final polishing.

Remineralizing agent

Contains fluoride to remineralize enamel after polishing, thus supporting caries prevention.

Integrated abrasion variability

ProphyCare PRO's high cleaning power converts into gentle polishing as the particles fragment.

ProphyCare PRO	
Starter	●
Very discoloured surfaces	
Thick plaque	
Discoloured teeth	●
Universal paste	●
Small amounts of plaque	●
Root surfaces	●
Slightly discoloured teeth	●
Fillings	●
Pediatric	●
Cervical toothbrush abrasion lesions	●
Implants	●
Final polishing	●
After air polishing treatment	●
To be used for polishing after scaling to relieve hypersensitivity after treatment	●

ProphyCare® PRO PRO, Lilac Universal

ProphyCare PRO is the one-step prophy paste which achieves outstanding cleaning and polishing results, as well as reduced treatment time. The same paste is used for stain removal and final polishing. ProphyCare PRO's high cleaning power converts into gentle polishing as the particles are fragmented. ProphyCare PRO offers effective stain removal properties at the beginning of treatment and decreases progressively to a lower RDA, giving the paste its final polishing properties at the end of treatment. ProphyCare's mild mint taste and appealing texture is pleasant for patients of all ages. It does not splatter, liquify or separate and is easy to rinse-off. It contains fluoride to remineralize enamel after polishing, thus supporting caries prevention. Available in single doses and tubes.

Prod No	Product Description
741111	ProphyCare PRO Tubes, 84g 60 ml
690125	ProphyCare PRO, Single Dose, 144 pcs, 2g



Interior webbing

Inner structure design minimizes splatter and increases friction for efficient stain removal.

Exterior ridges

Outer ridges optimize interproximal cleaning and splatter control. Soft flare improves subgingival access.

Soft edge

Flexible flare for improved subgingival access and efficient plaque removal.

Protective base

Wider base designed to prevent prophylactic paste from reaching the handpiece.



ProphyCare® Prophy Cup

Optimal cleaning, superior paste holding

The soft and flexible ProphyCare Prophy Cup flares and flexes to the contours of the teeth, aiding in selective stain removal and interproximal cleaning with less pressure and operator fatigue. Its size is optimized for easy access in hard-to-reach areas, without compromising the paste holding capacity. Its unique design is suitable for treatment on adults and children, and it can be used as a substitute for pointed polishers.

Prod No	Product Description
690131	ProphyCare Prophy Cups, screw-in, 144 pcs
690132	ProphyCare Prophy Cups, latch, 144 pcs






Rondell Blue and Rondell Red disclosing pellets

A convenient and cleaner way of applying disclosing solution. Simply rub the tooth surface with pre-soaked pellets and plaque will appear. Rondell blue is a bitonal plaque detector: new plaque, younger than 3 days, is coloured red, whilst older, more dangerous plaque is coloured blue. Both are easily brushed off, leaving a mild taste of peppermint. The dual-revealing colours are designed to educate patients and increase their compliance.

Prod No	Product Description
733552	Rondell, red, 100 pcs
733553	Rondell, blue, 100 pcs

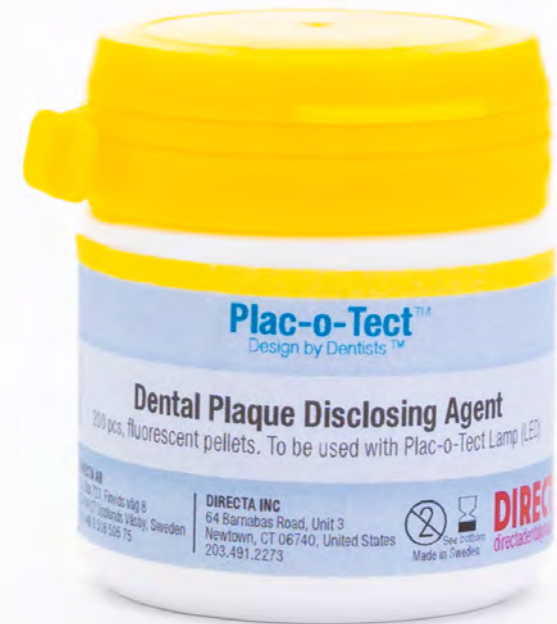

Plac-o-Tect™ Colors plaque made visible with UV-light

Plac-o-Tect Plaque disclosing agent is specially formulated so that it is only visible when exposed to UV-light. Plac-o-Tect is a simple and innovative way of demonstrating oral hygiene. The patients clearly see the results after the teeth have been cleaned. Invisible when the patient leaves the practice.

Prod No	Product Description
611450	Plac-o-Tect, refill, 1x 200 pc



Design By Dentists™

Directa's products are developed, tested and evaluated by our selected team of highly-qualified and renowned Swedish and international dental professionals. The aim is to make life easier for all dental professionals in their daily clinical work.

DIRECTA

DIRECTA AB

Finvidsväg 8,
194 47 Upplands Väsby, SWEDEN
+46 8 506 505 75

DIRECTA INC

64 Barnabas Road, Unit 3
Newtown, CT 06470, United States
203.491.2273

DIRECTA GERMANY

Hammergutstr. 11
92245 Küssmersbruck, GERMANY

info@directadental.com
www.directadental.com
www.directadentalgroup.com
www.facebook.com/directadental
www.instagram.com/directadental
www.youtube.com/user/DirectaAB

2022-03-14_REV01

DirectaDentalGroup

DIRECTA **topdental** **ORSING** **parkell** **physics**
forceps
RÖNVIG **TrollDental** **Sendoline** **CONTACTEZ**

