



OPTIMA EndoActivE+

Ultrasonic Endo Activator
Device BAE320



Ergonomic & comfortable design

Three Power Levels:

E1, E2, and E3, ranging from low to high power.

Four Time Modes:

Allows the ultrasonic to run for 10 seconds, 15 seconds, 20 seconds, or 10 minutes.

10s

15s

20s

10 min


Two Modes:

- 1) Stand-by ready to use mode
- 2) Irrigation Mode.

Highly functional design More efficient irrigation

Smart microcomputer control

EndoActiveE+ emits high-frequency ultrasonic vibrations that act on the handpiece head, making the tip vibrate. This vibration safely and effectively enhances the endo irrigation used in endodontics.

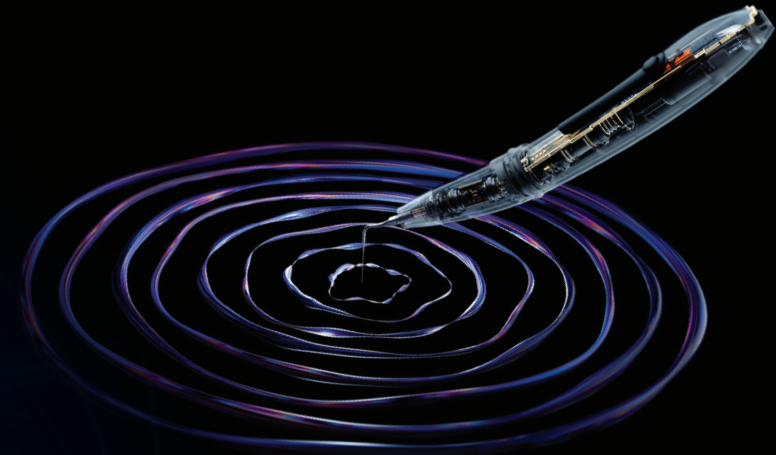

Wireless Control

User-friendly operation

Wireless Use & Wireless Charging Base: no plugs or cords.

Magnetic suction contact connects the handpiece to the base.

Compact and light for easy irrigation.



Clear LCD screen

Easier observation

Clear, large, easy to read display allows for better observation of the unit's power mode & timing.

Ring Switch

Control button

Turns the unit on and off and allows the user to switch between power levels and time modes along with the main power button.

Facilitates the endodontic treatment by allowing access to different angles and different tooth positions.



Stable performance

Appropriate vibration amplitude ensures consistent stable readings

High quality, high efficiency titanium alloy tips with low vibration ensure a reliable irrigation.



Intelligent tracking

Fully automatic frequency tracking system
Automatically search for best working state



Compatible tips



BAE70



BAE71



BAE72



BAE96



BAE99

High configuration: EN- BAE96×2, EN-BAE99×2, EN-BAE70×2, EN-BAE71×10, EN-BAE72×10

Compatible tips:EN- BAE70, EN-BAE71, EN-BAE72, EN-BAE96,EN- BAE99

For canal irrigation after the root canal preparation

EN-BAE70: Endochuck, used for install BAE71, BAE72, for root canal cleaning or irrigation
EN-BAE71, EN-BAE72: Polymer material without cutting ability

EN-BAE96: 20#, taper 2%, tip end length 21mm, Stainless steel

EN-BAE99: 25#, taper 2%, tip end length 18.5mm, Titanium alloy



 HENRY SCHEIN®

Henry Schein Australia

www.henryschein.com.au

Tel: 1300 65 88 22

Henry Schein NZ

www.henryschein.co.nz

0800 808 855