


CORDLESS ENDO MOTOR WITH CONTROL OF MOVEMENT, SPEED & TORQUE

GREAT SHAPING EXPERIENCE

- > Well-balanced in hand
- > Excellent visibility
- > Simple to operate
- > Adaptive, ideal for both left and right-handed practitioners

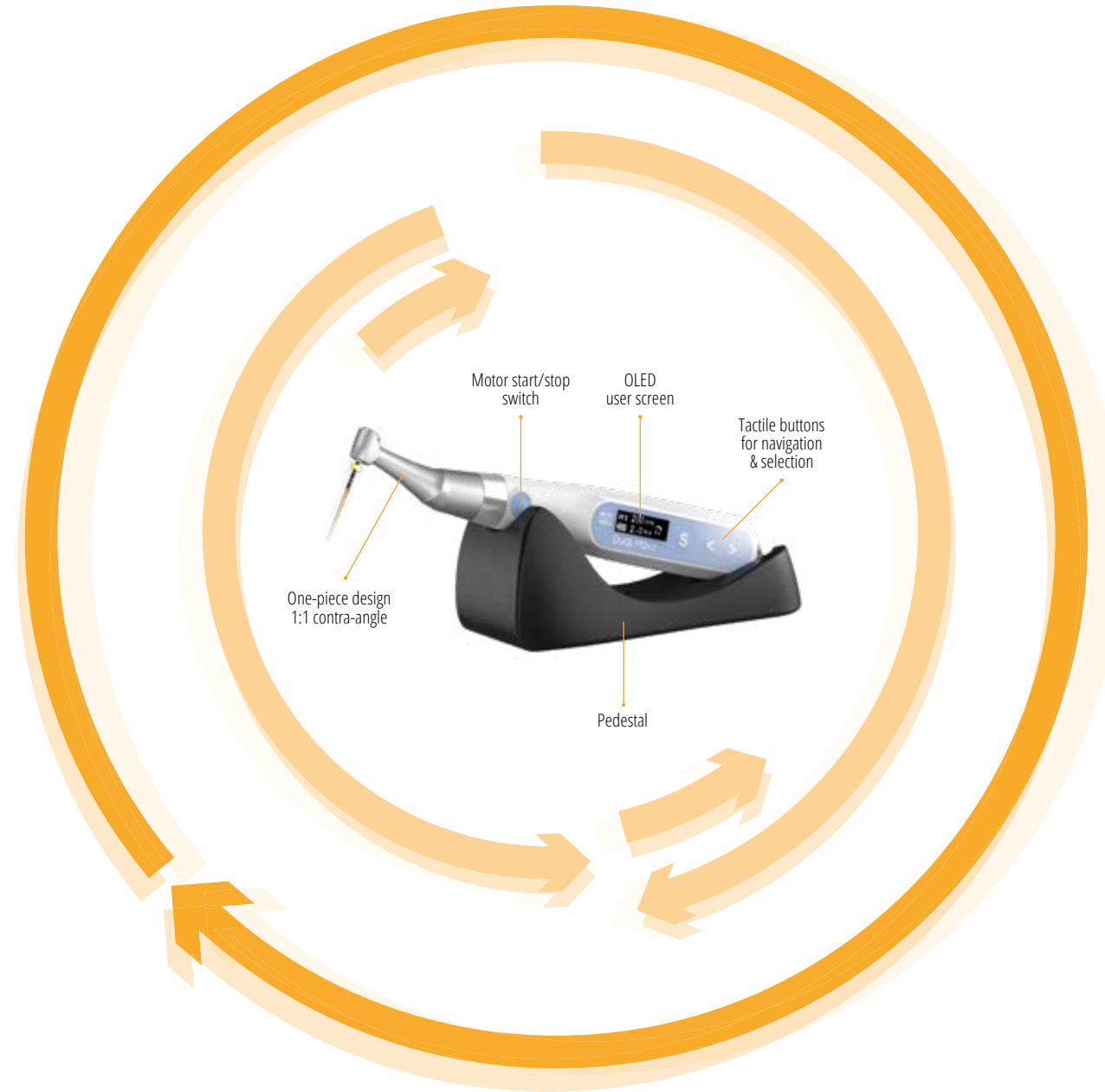

SAFE ENDO TREATMENT

- > Prevention of file separation
- > Speed and torque control

ADAPTED TO INDIVIDUAL NEEDS

- > Flexible settings with large number of parameters
- > Suitable for all NiTi files
- > Versatile, stand-alone use or combined with Dual Pex
- > Dual movement : continuous / reciprocating*

* Not available in some countries.

Dual·Move

GREAT SHAPING EXPERIENCE

> FREEDOM WHEN SHAPING!

- Slim and compact design
- Light and well-balanced handpiece



ADAPTED TO INDIVIDUAL NEEDS

> ALL LEVELS OF PRACTICE, ALL NiTi FILES!

- 6 customisable memory programmes with speed, rotation direction, torque limit
- Continuous rotation
- Reciprocating movement*





> TURN THE CONTRA-ANGLE TO THE PERFECT POSITION AND GET EXCELLENT VISIBILITY!

- 340° rotative contra-angle
- Miniature head and thin neck




Brand X



Dual Move

SAFE ENDO TREATMENT

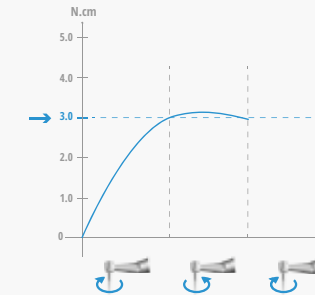

> PREVENT FILE SEPARATION!

- Auto-reverse on torque limit
- Alarms with adjustable volume level

> CONTROL YOUR PREPARATION WITH VISION AND SOUND!

- Ambidextrous display (rotation of 180°)
- Automatic locking of the file on the contra-angle
- Perfect position of start/stop switch for single hand operation
- User interface with short learning curve



> CONTROL YOUR TREATMENT!

- Reliable speed and torque control
- On-demand calibration function
- Instantaneous torque value on the screen and audible sounds



* Not available in some countries.