



EDGEAPEX^{HP}™ EDGEFLOW™ EDGEPACK™ EDGEAPEX™

EDGEPOWER

MORE POWER FOR ENDODONTICS

The Edge Power portfolio was designed to work seamlessly with our **NiTi and Obturation solutions.**





EDGEAPEX HP™

EDGEAPEX HP™ Designed to enhance precision and efficiency in endodontic procedures. Cordless and with an Integrated Apex Locator for maximum convenience.

Benefits

- Compatible with the majority of the existing file systems for better adaptability.
- Exceptional versatility with two movement options, continuous rotation and reciprocation.
- Cordless portable motor with wireless charging, auto on/off, and an integrated Apex Locator for maximum convenience.
- High-performance brushless motor with quiet operation and an extended service life.
- It adopts real-time feedback technology and dynamic torque contra angle, effectively preventing file separation.

Specifications:

- Great flexibility with a 360° rotation contra angle.
- Ultra-small head and slim neck for enhanced visibility and improved access to difficult areas.
- 10 available memories and pre-program settings for EdgeEndo file systems to save time.
- 4 available motions: CW / CCW / REC / ATR
- Multiple Apical functions
- Right-handed and left-handed adjustable display
- RPM range: 100 - 2500
- Torque range: 0.4 - 5.0 Ncm

The box includes

A motor handpiece, a charging base, a contra angle, a measuring wire, lip hook (x2), file clip (x4), touch probe (x2), power adapter, protective silicon cover (x2), spray nozzle, o-rings (x2), disposable sleeves (1 pack).

Item Code: **EN-EAPEXHP**



Motor Hand piece:
EN-EAPEXRU
Charging Base: **EN-EAPEXCB**

Main accessories

Product	Item Code	Image
MEASURING WIRE	EN-EAPEXNMW	
TOUCH PROBE	EN-EAPEXTP	
LIP HOOK	EN-EAPEXLH	
CONTRA ANGLE HANDPIECE CA 161	EN-EAPEXCA	
FILE CLIP	EN-EAPEXFC	

Warranty: 2 years for handpiece motor, charging base, and power adapter. 1 year for contra angle attachment. 6 months for accessories.



EDGEFLOW™

360° Rotating with excellent gutta-percha fluidity for a uniform, and smooth, bubble-free output ensuring seamless fillings.

Benefits

- Excellent gutta-percha fluidity for a uniform, bubble-free output ensuring seamless fillings.
- Wireless and electric design for maximum comfort and ease of use.
- Dual battery system provides extended endurance, allowing for approximately 100 root canal fillings.
- Precise and intuitive display to better monitor the remaining gutta-percha.

Specifications:

- Auto heating stop after timeout; Auto shutdown after timeout
- Charging base with battery
- Temperature monitoring design
- Reusable injection needles available in various sizes and lengths
- Preset temperatures: 100°C, 120°C, 150°C, 180°C, 200°C
- Injection Speeds: Slow, Medium, Fast



Main Unit: **EN-EEPTORU**
Charging Base: **EN-EEPTOCH**

The box includes

A main unit, Charging Stand, Power Adapter, Battery (x2), Heat Insulation Hood, Wrench of Needle, Gutta-Percha Ring, Injection Needle (x4) E23G-M (1) - E23G-L (1) - E25G-L (2).

Item Code: **EN-EEPTO**

Main Accessories

Product	Item Code	Image
WRENCH	EN-EEPTOW	
LITHIUM-ION BATTERY	EN-EEPTOB	
THERMAL PROTECTOR CAP	EN-EEPTOHIH3	

Tips

Product	Item Code	Image
*E20G-S INJECTION NEEDLE	EN-EE20GSIN	
*E20G-M INJECTION NEEDLE	EN-EE20GMIN	
*E20G-L INJECTION NEEDLE	EN-EE20GLIN	
E23G-M INJECTION NEEDLE	EN-EE23GMIN	
E23G-L INJECTION NEEDLE	EN-EE23GLIN	
*E25G-M INJECTION NEEDLE	EN-EE25GMIN	
E25G-L INJECTION NEEDLE	EN-EE25GLIN	

Warranty: 1 year for device.



EDGEPACK™

EDGEPACK™ Features high-precision temperature control with pre-set working temperatures, ensuring optimal compatibility with various gutta-percha types for precise and reliable obturation.

Benefits



- Wireless Gutta-Percha Obturation Device for seamless, cordless operation.
- Ergonomic and ambidextrous design for enhanced comfort and usability.
- 0.35mm pre-bend tip for improved access and better penetration into narrow root canals.
- High-precision temperature control for consistent and reliable obturation.
- Efficient GP cutting system.
- Compact size for easy cleaning and convenient storage.
- Visual and auditory prompts for enhanced safety during procedures.

Specifications:

- Preset temperature options: 150°C, 180°C, 200°C, and 230°C.
- Advanced time-related features:
- Ultra-fast warm-up time of 0.2 seconds.
- Auto shut-off after 10 seconds for energy efficiency.
- Long-lasting battery life—provides up to 4 hours per charge.


The box includes

A Charging Stand, Power Adapter, WP4004 Tip, WP4504 Tip, WP5506 Tip, WP5508 Tip, Heat Insulation Sleeve of Tips (x4).









Main Unit: **EN-EEDPURU**

Item Code: **EN-EEDPU**

Main Accessories

Product	Item Code	Image
CHARGING BASE	EN-EEDPUCS	

Tips

Product	Item Code	Image	Product	Item Code	Image
WP4004 WORK TIP	EN-EEDPUT4004		EDGEPACK WP 5508 WORK TIP	EN-EEDPUT5508	
WP4504 WORK TIP	EN-EEDPUT4504		*EDGEPACK WP 5508L WORK TIP	EN-EEDPUT5508L	
*WP5004 WORK TIP	EN-EEDPUT5004		*EDGEPACK WP 5510 WORK TIP	EN-EEDPUT5510	
WP5506 WORK TIP	EN-EEDPUT5506		*EDGEPACK WP 6012 WORK TIP	EN-EEDPUT6012	

Warranty: 1 year for device and tips.



EDGEAPEX™

EDGEAPEX™ Effectively locates the apical constriction, enabling precise determination of root canal working length for optimal treatment outcomes.



Main Unit: **EN-EEALURU**
Charging Base: **EN-EEALUB**

Benefits

The EdgeApex is an apex locator and a pulp vitality tester.






It utilises cutting-edge multi-frequency technology combined with a highly efficient digital signal processor, delivering optimal accuracy and featuring a high resolution large screen.

Pulp vitality tests diagnose pulpitis, determine the need for root canal treatment or other dental procedures, and assess overall tooth health.

The box includes

- Charging Base (1)
 - Measuring Wire (1)
 - File Clip (4)
 - Lip Hook (5)
 - Touch Probe (2)
 - Pulp Tester Probe (2)
 - Tester (1)
 - Lithium-ion Battery (1)
 - Power Adapter (1)
- Enlarged display of apical foramen (1)
- Item Code: **EN-EEALU**

Main Accessories

Product	Item Code	Image
MEASURING WIRE	EN-EEALUMW	
TOUCH PROBE	EN-EEALUMF	
LIP HOOK	EN-EEALULH5	
ENDO RULER	EN-EEALUER	
FILE CLIP	EN-EEALUFC5	
PULP TESTER PROBE	EN-EEALUPTP	

Warranty: 1 year for device. 6 months for accessories.

EDGEPOWER WORKFLOW

Explore our extensive range of endodontics products on our website www.edgeendo.com.au or www.henryschein.co.nz/edgeendo

1. ACCESS

- EDGEAPEX HP™
- EDGE R GLIDE
- EDGE APEX™ Locator



2. SHAPING

- EDGEAPEX HP™
- EDGEOne R




3. IRRIGATION

- Paper Points



4. OBTURATION

- EDGEBIOCERAMIC™ Sealer
- Gutta Percha Points
- EDGEFlow™
- EDGEPack™

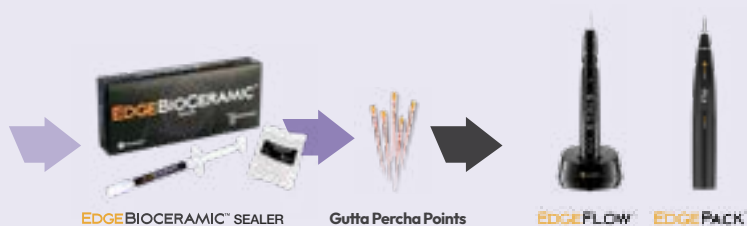

 HENRY SCHEIN®

AUSTRALIA

1300 65 88 22 | Customer.Care@henryschein.com.au | www.henryschein.com.au

NEW ZEALAND

0800 808 855 | info@henryschein.co.nz | www.henryschein.co.nz