

SciCan HYDRIM C61W G4

Instrument Washer



Automated Instrument Washing the first step in the sterilization process



Efficient

The time saved by automated washing compared to manual cleaning is approximately 1 hour of labour for every nine procedural set-ups¹. HYDRIM C61W G4 cycle times are as fast as 40 minutes from start to dry. HYDRIM ensures that every set of instruments receives the same validated cleaning process.



G4 Technology

There is nothing like SciCan's G4 Technology. The large touch screen gives in-depth cycle information and allows users to view current operations.

From the networking portal, users can access cycle history and validation information. If there is ever a cycle fault, the G4 software will send an email to the office and the service technician, allowing remote troubleshooting.



Safe

Avoid the many steps involved with manual washing or using ultrasonic cleaners. Steps such as pre-soaking, scrubbing and drying instruments before sterilization increase the exposure to blood and body fluids and the risk of sharps injuries.

Automated instrument washing reduces physical handling of contaminated instruments. Using the HYDRIM G4 instrument washer is the safer and more effective option, decreasing the risk of exposure and injuries. HYDRIM G4 saves time and money—making it a worthy investment for today's standard of care.



Capacity

HYDRIM C61W G4 has a 6 full-size cassette capacity, allowing offices to turnover 120 instruments in as little as 40 minutes, start to dry.

¹ Tholen, Mark. "Increase production and decrease expenses in sterilization." Dental Economics, Volume 99, Issue 10, October 1, 2009.

HYDRIM C61W G4 Specifications

Unit Size

Dimensions (W x H x D):	598 x 520 x 526mm / 23.5" x 20.5" x 20.7"
Depth with Door Open:	829 mm / 32.6"
Weight:	44 kg / 97 lb

Cycle Times*

Rinse and Hold	5 min
Regular Wash + Dry	26 min + 10 min
Heavy Duty Wash + Dry	30 min + 10 min

Water Consumption Per Cycle

Rinse and Hold	7.6 L / 1.7 Gal
Regular Wash	11.5 L / 3 Gal
Heavy Duty Wash	11.5 L / 3 Gal

Electrical Connection

North America	208-240 V, 12 A, 60 Hz
Running Noise	65 dB

Water and Drain Connections

Hot and Cold Water Tap Connection	3/4" Garden Hose, Male
Recommended Hot Water Temperature	60° C / 140° F
Cold Water Temperature	< 30° C / 86° F
Inlet Water Pressure	1-10 bar / 14.5-145 psi
Drain	2 cm / 3/4"

*Actual cycle times will vary with incoming water temperature and pressure.

Recommended Accessories

Part Number	Description
01-113251	6-Cassette Rack
01-113252	4 Cassette Rack
01-113253	3-Cassette, 2-Basket Rack
01-113254	10-Tray, 2-Basket Rack
01-113257	Trolley Center Support
01-113545	Large Basket
01-113546	Hygiene Basket
01-113547	Basket with Hinged Lid
CS-HIPC-U	HIP ULTRA Cleaning Solution



VistaPure Water Filtration

Provides two grades of high-quality water for autoclaves, dental bottles, ultrasonic cleaners and the rinse cycle of instrument washers.

V3000

COLTENE International Dental Group

Manufactured for:
Dent4You AG
 Bahnhofstrasse 2
 CH-9435 Heerbrugg

Distributed by:
SciCan Ltd.
 1440 Don Mills Rd.
 Toronto, ON
 M3B 3P9 Canada
 T +1 416-445-1600
 customerservice@scican.com


www.scican.com

HYDRIM VistaPure Connection

Protect expensive instruments and equipment by using the VistaPure system's hyperfiltered water for the final-rinse cycle of the HYDRIM for a spot-free finish.

V3010

VistaPure is manufactured by:
Crosstex International, Inc.
 1733 W. Parkside Lane
 Phoenix, AZ 85027 / USA
 T +1 800-531-2221
 VistaResearchGroup.com

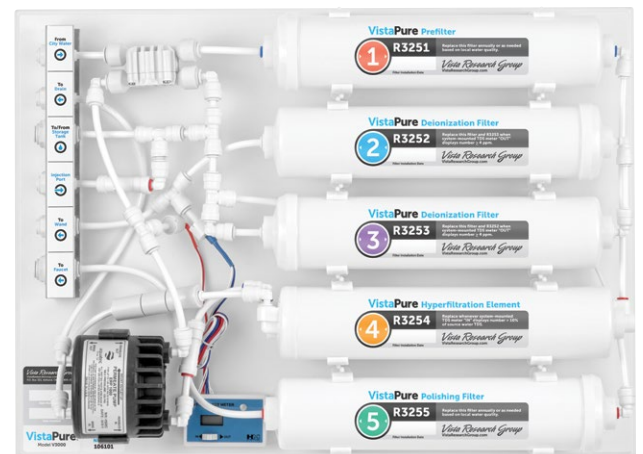

VistaPure Water Filtration shown above.