



OSTEOGENICS

BIOMEDICAL | EST. 1996



Regeneration Products Catalog

osteogenics.com

Welcome to Osteogenics Biomedical

We are Osteogenics Biomedical, makers of Cytoplast™. Established in 1996 with a goal to create a more predictable alternative to Gore-Tex® membranes, we have grown to be a leader in barrier membrane and PTFE suture technologies in the United States. After 20 years of product development focused on surgical predictability, we are expanding globally. We encourage you to try Cytoplast™ regenerative products to see why thousands of surgeons rely on us. We guarantee your satisfaction – or your money back. *To find the distributor nearest you, go to www.osteogenics.com/GlobalNetwork.*

Table of Contents

4 **Zcore™ Porcine Xenograft Particulate**

Resorbable GTR Barrier Membranes

6 Cytoplast™ RTM Collagen

7 Zmatrix™ Porcine Peritoneum Collagen Membrane

8 Vitala® Porcine Derived Collagen

Non-Resorbable GTR Barrier Membranes

9 Cytoplast™ TXT-200 & TXT-200 Singles

10 Cytoplast™ Ti-250 & Ti-150 Titanium-Reinforced

14 • RPM™ Reinforced PTFE Mesh

18 • **Cytoplast™ PTFE Suture**

Pro-Fix™ – Precision Fixation System

19 Pro-Fix™ Membrane Fixation


19 Pro-Fix™ Individual Components

20 Pro-Fix™ Tenting

21 Pro-Fix™ Bone Fixation

22 **Selection of Applicable References**

• **New Items Available**

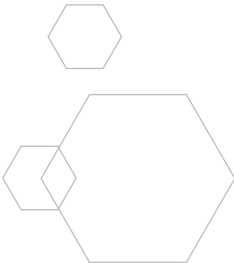



Zcore™ Porcine Xenograft Particulate

.25 mm - 1.0 mm Particle Size

ZS050	0.5 cc
ZS100	1.0 cc
ZS200	2.0 cc
ZS400	4.0 cc



Zcore™ Porcine Xenograft Particulate

1.0 mm - 2.0 mm Particle Size

ZL100	1.0 cc
ZL200	2.0 cc



Zcore™ Porcine Xenograft Particulate in Syringe

.25 mm - 1.0 mm Particle Size

ZY025	0.25 cc
ZY050	0.5 cc



Features & Benefits of Zcore™

Zcore™ is an osteoconductive, porous, anorganic bone mineral with a carbonate apatite structure derived from porcine cancellous bone.

Interconnecting pores

Interconnecting macroscopic and microscopic porous structure supports the formation and ingrowth of new bone

88% to 95% void space

88% to 95% Void Space: hyper-porosity of porcine cancellous matrix and intra-particle space facilitated by rough particle morphology reduce bulk density of the graft, allowing greater empty space for new bone growth*

Porcine cancellous bone


Derived from porcine cancellous bone, eliminating risk of BSE transmission

Processed using minimal heat

Heat treated to an optimal temperature that ensures a degree of crystallinity¹ consistent with native bone mineral to allow for remodeling of the healing bone


*0.25 mm – 1.0 mm particle size = 88% void space, 1.0 mm – 2.0 mm = 95% void space

1. Li ST, Chen HC, Yuen D. Isolation and Characterization of a Porous Carbonate Apatite From Porcine Cancellous Bone. Science, Technology, Innovation, Aug. 2014: 1–13.

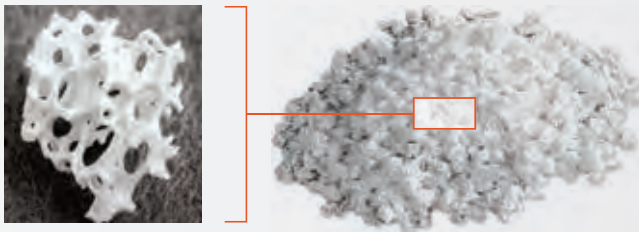


SEM of Processed Human Bone

Magnification x50

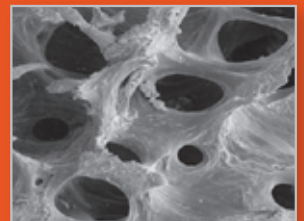

SEM of Zcore™ Porcine Xenograft Particulate


Magnification x50

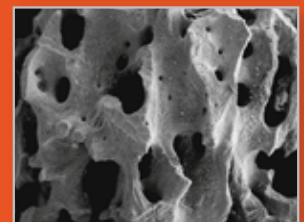

Cytoplast™ RTM Collagen

Type I bovine collagen membrane

shown actual size.

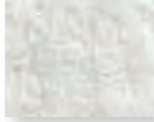

15 mm x 20 mm

RTM1520 (2 membranes per box)



20 mm x 30 mm

RTM2030 (2 membranes per box)



30 mm x 40 mm

RTM3040 (2 membranes per box)

Features & Benefits

Manufactured from highly purified type I bovine achilles tendon

Safe for the patient

26 – 38 week resorption time

Long predictable resorption time limits the risk of particle loss due to premature resorption

High tensile strength


You can suture or tack the membrane in place without tearing

Cell occlusive

Prevents epithelial down growth

Optimized flexibility

Stiff enough for easy placement, yet easily drapes over ridge

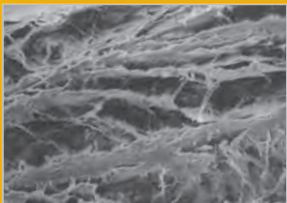

Multi-layer construction allows tissue integration into outer layer, while preventing direct passage of bacteria and epithelial cells.

“...I am impressed with its *handling*, but most importantly, I am impressed with its *results*.”

Jerald Rosenberg, DMD; Periodontist

Zmatrix™

Porcine peritoneum collagen membrane



A perfectly soft consistency that drapes without the usual self-adherence experienced with other natural collagen membranes

shown actual size.

15 mm x 20 mm

ZM1520



20 mm x 30 mm

ZM2030



30 mm x 40 mm

ZM3040



Features

Extracellular Components

Processed to preserve extracellular components including laminin, fibronectin, elastin, and glycosaminoglycans*

Easy to Handle

Designed to drape without adhering to itself

Elastic

Natural peritoneum collagen structure allows for elasticity

Natural, Native Collagen Membrane

*Zmatrix™ is a natural, native collagen membrane; cross-linking chemicals and agents are unnecessary. Proprietary processing technology allows preservation of collagen as well as extracellular components including laminin, fibronectin, elastin, and glycosaminoglycans.**

*Hoganson DM, Owens GE, O'Doherty EM, Bowley CM, Goldman SM, Harilal DO, Neville CM, Kronengold RT, Vacanti JP. Preserved extracellular matrix components and retained biological activity in decellularized porcine mesothelium. Biomaterials. 2010, 27: 6934–6940.



Vitala®

Porcine pericardium collagen membrane | Substantially resorbed in 26 weeks

shown actual size.



10 mm x 10 mm
VIT1010

vitala mini



15 mm x 20 mm
VIT1520




13 mm x 25 mm
VIT1325



20 mm x 30 mm
VIT2030



30 mm x 40 mm
VIT3040

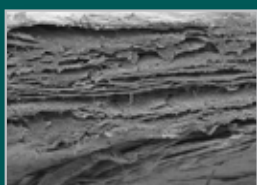

1000x magnification



Excellent tensile strength



Supple and flexible

Features & Benefits

Natural

Manufactured using a proprietary protocol designed to maintain the natural, microporous, 3-layered architecture of the tissue without the need for cross-linking chemicals and agents

Durable

Designed to resist tearing during placement, Vitala® is naturally strong

Adaptable

The natural collagen structure provides a unique combination of supple handling and ideal defect adaptability. Because both sides are smooth, either side may be placed against the defect

Cytoplast™ TXT-200 & TXT-200 Singles

Micro-textured, high-density PTFE membrane

Most popular
membrane for
socket grafting

TXT-200 Singles 12 mm x 24 mm

TXT1224-1 (1 membrane per box)

TXT1224 (10 membranes per box)

shown actual size.



TXT-200 25 mm x 30 mm

TXT2530-1 (1 membrane per box)

TXT2530 (4 membranes per box)



Features & Benefits

Non-Resorbable

Won't resorb prematurely – you dictate healing time

100% Dense (non-expanded) PTFE

Impervious to bacteria (pore size less than 0.3 µm)

Data on file

Purposely leave the membrane exposed

Preservation of the soft tissue architecture and
keratinized mucosa


Soft tissue attaches, but doesn't grow through the membrane

Exposed membrane allows for non-surgical removal;
no anesthesia required

Hexagonal dimples increase surface area

Designed to increase membrane stabilization

The patented Regentex™ surface helps stabilize the membrane and the soft tissue flap. Hexagonal surface dimples provide a textured surface that increases the area available for cellular attachment without increasing porosity.
U.S. Patent # 5,957,690

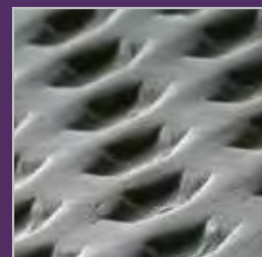

"I always know, *in advance*, the results of my bone grafting when I use Cytoplast™ TXT-200 as a membrane. *Why bother with other membranes?*"

Mark Cohen, DDS; Periodontist




Cytoplast™ Titanium-Reinforced

Titanium-reinforced, high-density PTFE membrane

	Ti-250 (250 µm thick)	Ti-150 (150 µm thick)		Versatile Rectangular Shapes
				<p>These configurations can be trimmed to fit a variety of defects.</p> <p>Shown actual size.</p>
ANL 12 mm x 24 mm Designed for narrow single-tooth extraction sites, especially where one bony wall is missing	Ti250ANL-1 Ti250ANL-2	Ti150ANL-1 Ti150ANL-2	(1 membrane per box) (2 membranes per box)	
ANL30 12 mm x 30 mm Designed for narrow single-tooth extraction sites, especially where one bony wall is missing	Ti250ANL30-1 Ti250ANL30-2		(1 membrane per box) (2 membranes per box)	
PS 20 mm x 25 mm Designed for large extraction sites and limited ridge augmentation	Ti250PS-1 Ti250PS-2	Ti150PS-1 Ti150PS-2	(1 membrane per box) (2 membranes per box)	
PL 25 mm x 30 mm Designed for large bony defects, including ridge augmentation	Ti250PL-1 Ti250PL-2	Ti150PL-1 Ti150PL-2	(1 membrane per box) (2 membranes per box)	




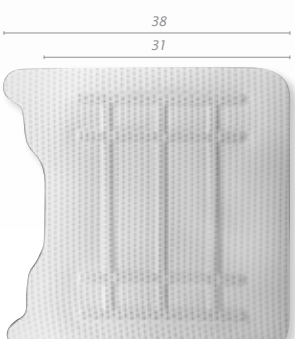


**Ti-150 membranes are 40% thinner than Ti-250 membranes, providing clinicians another handling option in Cytoplast™ Titanium-Reinforced Membranes.*

	Ti-250 (250 µm thick)	Ti-150 (150 µm thick)		Versatile Rectangular Shapes
				<p>These configurations can be trimmed to fit a variety of defects.</p> <p>Shown actual size.</p>
XL 30 mm x 40 mm Designed for very large bony defects, including ridge augmentation	Ti250XL-1	Ti150XL-1	(1 membrane per box)	
	Ti250XL-2	Ti150XL-2	(2 membranes per box)	
XLK 30 mm x 40 mm Designed for very large bony defects, including ridge augmentation	Ti250XLK-1	Ti150XLK-1	(1 membrane per box)	
	Ti250XLK-2	Ti150XLK-2	(2 membranes per box)	
K2 40 mm x 50 mm Designed for the largest bony defects, including ridge augmentation	Ti250K2-1	Ti150K2-1	(1 membrane per box)	
	Ti250K2-2	Ti150K2-2	(2 membranes per box)	

Cytoplast™ Titanium-Reinforced

Titanium-reinforced, high-density PTFE membrane

	Ti-250 (250 µm thick)	Ti-150 (150 µm thick)		Interproximal Shapes
				<p>These configurations are designed to fit between existing teeth.</p> <p>Dimensional measurements shown in mm Width measurements noted at widest point and narrowest point. Shown actual size.</p>
AS 14 mm x 24 mm Designed for single-tooth extraction sites, especially where one or more bony walls are missing	Ti250AS-1 Ti250AS-2	Ti150AS-1 Ti150AS-2	(1 membranes per box) (2 membranes per box)	
ATC 24 mm x 38 mm Designed for large extraction sites, including ridge augmentation	Ti250ATC-1 Ti250ATC-2	Ti150ATC-1 Ti150ATC-2	(1 membranes per box) (2 membranes per box)	
PTC 38 mm x 38 mm Designed for large bony defects, including ridge augmentation	Ti250PTC-1 Ti250PTC-2	Ti150PTC-1 Ti150PTC-2	(1 membranes per box) (2 membranes per box)	
PD 38 mm x 38 mm Designed for large bony defects, including distal extension of the posterior ridge	Ti250PD-1 Ti250PD-2	Ti150PD-1 Ti150PD-2	(1 membranes per box) (2 membranes per box)	

Cytoplast™ Titanium-Reinforced

Titanium-reinforced, high-density PTFE membrane

	Ti-250 (250 µm thick)	Ti-150 (150 µm thick)		Shapes with Fixation Points
BL 17 mm x 25 mm Designed for large buccal defects	Ti250BL-1	Ti150BL-1	(1 membranes per box)	<p>These configurations are designed with fixation points outside of the defect area.</p> <p>Dimensional measurements shown in mm Width measurements noted at widest point and narrowest point. Shown actual size.</p>   
	Ti250BL-2	Ti150BL-2	(2 membranes per box)	
PST 36 mm x 25 mm Designed for large extraction sites and limited ridge augmentation in the anterior maxilla	Ti250PST-1	Ti150PST-1	(1 membranes per box)	
	Ti250PST-2	Ti150PST-2	(2 membranes per box)	
PLT 41 mm x 30 mm Designed for large bony defects, including ridge augmentation in the anterior maxilla	Ti250PLT-1	Ti150PLT-1	(1 membranes per box)	
	Ti250PLT-2	Ti150PLT-2	(2 membranes per box)	
	Ti-250 (250 µm thick)	Ti-150 (150 µm thick)		Perio Shapes
AP 13 mm x 19 mm Designed for periodontal defects in the anterior	Ti250AP-1	Ti150AP-1	(1 membranes per box)	<p>These configurations are designed for grafting perio defects.</p> <p>Shown actual size.</p> 
	Ti250AP-2	Ti150AP-2	(2 membranes per box)	
PP 13 mm x 18 mm Designed for periodontal defects in the posterior	Ti250PP-1	Ti150PP-1	(1 membranes per box)	
	Ti250PP-2	Ti150PP-2	(2 membranes per box)	



RPM™
Reinforced PTFE mesh

NEW

Versatile Rectangular Shapes

These configurations can be trimmed to fit a variety of defects.

Shown actual size.

- **PS**
20 mm x 25 mm
Designed for large extraction sites and limited ridge augmentation
- **PL**
25 mm x 30 mm
Designed for large bony defects, including ridge augmentation
- **XL**
30 mm x 40 mm
Designed for very large bony defects, including ridge augmentation
- **XLK**
30 mm x 40 mm
Designed for very large bony defects, including ridge augmentation
- **XLKM (mandible)**
30 mm x 40 mm
Designed for very large bony defects, including mandibular ridge augmentation *NOTE: Non-perforated region is designed for lingual aspect*

RPM250PS



RPM250PL



RPM250XL



RPM250XLK



RPM250XLKM





RPM's unique circular macroporous design allows for direct contact between the bone graft and periosteum, allowing naturally occurring revascularization and infiltration of cells into the bone graft.

- **K2**
40 mm x 50 mm
Designed for the largest bony defects, including ridge augmentation

RPM250K2

Versatile Rectangular Shapes





- **BL**
17 mm x 25 mm
Designed for large buccal defects

RPM250BL

Shapes with Fixation Points

These configurations are designed with fixation points outside of the defect area.

*Dimensional measurements shown in mm
Width measurements noted at widest point and narrowest point. Shown actual size.*

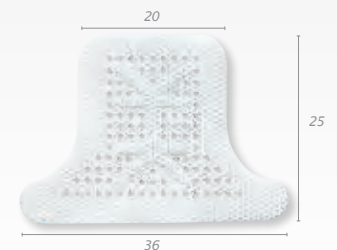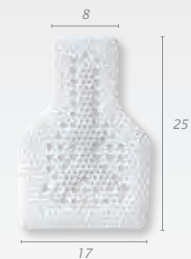

- **PST**
36 mm x 25 mm
Designed for large extraction sites and limited ridge augmentation in the anterior maxilla

RPM250PST

- **PLT**
41 mm x 30 mm
Designed for large bony defects, including ridge augmentation in the anterior maxilla

RPM250PLT

Interproximal Shapes

These configurations are designed to fit between existing teeth.

Dimensional measurements shown in mm
Width measurements noted at widest point and narrowest point. Shown actual size.

• **ATC**

24 mm x 38 mm

Designed for large extraction sites, including ridge augmentation

RPM250ATC

• **ATCM (mandible)**

24 mm x 38 mm

Designed for large extraction sites, including mandibular ridge augmentation
NOTE: Non-perforated region is designed for lingual aspect


RPM250ATCM

• **PTC**

38 mm x 38 mm

Designed for large bony defects, including ridge augmentation

RPM250PTC

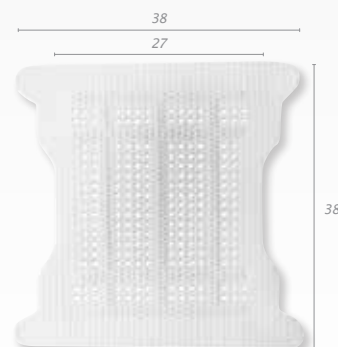• **PTCM (mandible)**

38 mm x 38 mm

Designed for large bony defects, including mandibular ridge augmentation
NOTE: Non-perforated region is designed for lingual aspect

RPM250PTCM



Interproximal Shapes


These configurations are designed to fit between existing teeth.

*Dimensional measurements shown in mm
Width measurements noted at widest point and narrowest point. Shown actual size.*

- **PD**
38 mm x 38 mm

Designed for large bony defects, including distal extension of the posterior ridge

RPM250PD

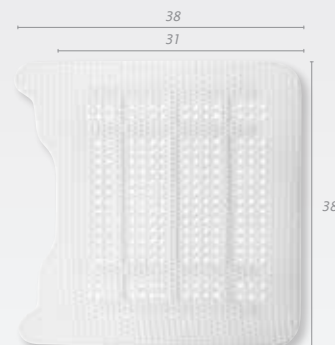

- **PDMR (mandible right)**
38 mm x 38 mm

Designed for large bony defects, including distal extension of the right posterior mandibular ridge
NOTE: Non-perforated region is designed for lingual aspect


RPM250PDMR



- **PDML (mandible left)**
38 mm x 38 mm

Designed for large bony defects, including distal extension of the left posterior mandibular ridge
NOTE: Non-perforated region is designed for lingual aspect

RPM250PDML

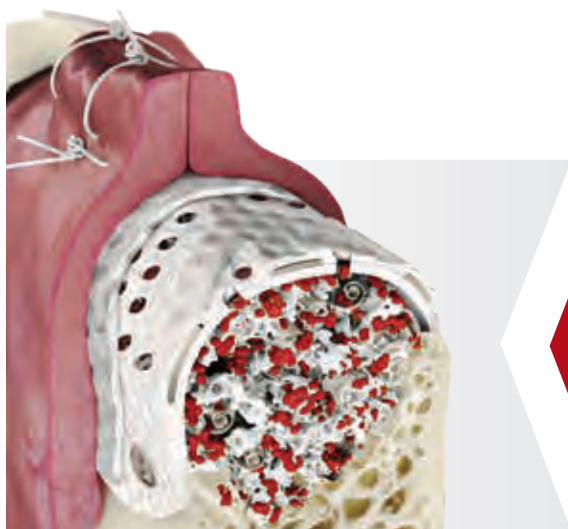

Circular Macropores allow direct contact between bone graft and periosteum, allowing naturally occurring revascularization and infiltration of cells into the bone graft

Titanium Frame maintains space essential for horizontal and vertical ridge augmentation

PTFE Mesh easily conforms to tissue contours

Cytoplast™ PTFE Suture

The soft monofilament suture

300 Series Stainless Steel Needles

All Cytoplast™ PTFE Sutures now have 300 series stainless steel needles, the gold standard material for suture needles. Tests comparing the new needles to previous needles show a substantial increase in needle strength, initial needle sharpness, and sustained needle sharpness. Tests show that the new 300 series needles are less likely to bend, require less force to penetrate, and maintain sharpness longer. Additionally, CS0618RC and CS06PREM now have longer (121% and 41%, respectively) and geometrically finer precision cutting edges. Data on file

Cytoplast™ undyed 19 mm precision RC 2/0 USP	CS0418
Cytoplast™ undyed 16 mm precision RC 3/0 USP	CS0518
Cytoplast™ undyed 19 mm precision RC 3/0 USP	CS051819

Cytoplast™ undyed 16 mm RC black needle 3/0 USP	CS0518BK
Cytoplast™ undyed 19 mm RC black needle 3/0 USP	CS051819BK

Cytoplast™ undyed 13 mm TP 4/0 USP	CS0618PERIO
Cytoplast™ undyed 13 mm precision RC 4/0 USP	CS0618PREM
Cytoplast™ undyed 16 mm precision RC 4/0 USP	CS0618RC
• Cytoplast™ undyed 13 mm precision RC 5/0 USP	CS071813
• Cytoplast™ undyed 16 mm precision RC 5/0 USP	CS071816

NEW

Needle Code Detail

RC 3/8 Circle Reverse Cutting



TP 1/2 Circle Round-Bodied



Features & Benefits

100% Medical Grade PTFE

Biologically inert

Monofilament

Doesn't wick bacteria

Soft (not stiff)

Comfortable for patients

Little to no package memory

Excellent handling, knots securely

Non-resorbable

Keeps the surgical site reliably closed

Pro-Fix™ Membrane Fixation

Precision Fixation System

Pro-fix™ Membrane Fixation Screws are designed as an attractive alternative to using tacks for membrane stabilization. Easy pick-up, solid stability of the screw during transfer to the surgical site, and easy placement make membrane fixation fast and easy.

Tray and organizer dial are designed to store all Pro-fix™ components including up to 100 membrane fixation, tenting, and bone fixation screws

Blades are designed to work universally with all Pro-fix™ membrane fixation, tenting, and bone fixation screws



Membrane Fixation Kit

PFMK20

- (1) Autoclavable Tecapro™ storage tray w/ screw organizer dial
- (1) Stainless steel driver handle
- (1) 76 mm cruciform driver blade
- (1) 56 mm cruciform driver blade
- (20) 1.5 x 3.0 mm self-drilling membrane fixation screws



Self-Drilling Membrane Fixation Screws

1.5 mm x 3.0 mm  actual size


- | | |
|-----------|---------|
| 5 screws | PFMF-5 |
| 10 screws | PFMF-10 |
| 20 screws | PFMF-20 |



Individual Components

- | | |
|---|--------|
| Stainless Steel Driver Handle | PFDH |
| 76 mm Cruciform Driver Blade | PFDB |
| 56 mm Cruciform Driver Blade | PFDB56 |
| Contra Angle Blade | PFDBCA |
| (24 mm long; 10 mm exposed distal length) | |
| Autoclavable Tecapro™ storage tray | PFT |
| 1.2 mm diam. Latch Type Pilot Drill | PFPD |



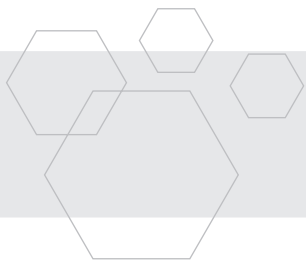

Pro-Fix™ Tenting

Precision Fixation System

Tenting Kit

PFTK12

- (1) Autoclavable Tecapro™ storage tray w/ screw organizer dial
 - (1) Stainless steel driver handle
 - (1) 76 mm cruciform driver blade
 - (1) 56 mm cruciform driver blade
 - (4) 1.5 x 3.0 mm self-drilling tenting screws (7 mm total length: see below)
 - (4) 1.5 x 4.0 mm self-drilling tenting screws (8 mm total length: see below)
 - (4) 1.5 x 5.0 mm self-drilling tenting screws (9 mm total length: see below)
- For individual Pro-Fix™ driver and container components, see page 19.

Pro-fix™ Tenting Screws are designed with a self-drilling tip, polished neck, and broader head to maintain space under resorbable and non-resorbable membranes in horizontal and vertical bone regeneration procedures.

Self-Drilling Tenting Screws

1.5 mm x 3.0 mm

3.0 mm polished neck + 4.0 mm threaded portion = 7 mm total length

- | | |
|----------|--------|
| 1 screw | PFT3 |
| 5 screws | PFT3-5 |

 actual size



1.5 mm x 4.0 mm

4.0 mm polished neck + 4.0 mm threaded portion = 8 mm total length

- | | |
|----------|--------|
| 1 screw | PFT4 |
| 5 screws | PFT4-5 |

 actual size

1.5 mm x 5.0 mm

5.0 mm polished neck + 4.0 mm threaded portion = 9 mm total length

- | | |
|----------|--------|
| 1 screw | PFT5 |
| 5 screws | PFT5-5 |

 actual size

Fully Threaded Tenting Screws

1.5 mm x 8.0 mm

- | | |
|---------|------|
| 1 screw | PFT8 |
|---------|------|

 actual size

1.5 mm x 10.0 mm


- | | |
|---------|-------|
| 1 screw | PFT10 |
|---------|-------|

 actual size



Pro-Fix™ Bone Fixation

Precision Fixation System

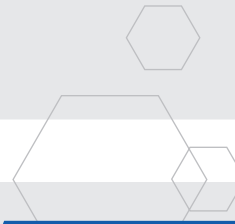

Bone Fixation Kit

PFBK12

- (1) Autoclavable Tecapro™ storage tray w/ screw organizer dial
- (1) Stainless steel driver handle
- (1) 76 mm cruciform driver blade
- (1) 56 mm cruciform driver blade
- (1) 1.2 mm diameter latch type pilot drill
- (2) 1.5 x 8 mm bone fixation screws
- (4) 1.5 x 10 mm bone fixation screws
- (4) 1.5 x 12 mm bone fixation screws
- (2) 1.5 x 14 mm bone fixation screws

For individual Pro-Fix™ driver and container components, see page 19.

Pro-fix™ Bone Fixation Screws are designed with finer pitched, self-tapping threads that give the screws greater clamping force while using less driver torque. The screws' threads are equipped with a cutting flute that allows for easier insertion into harder bone. The screws are placed into a 1.2 mm pre-drilled pilot hole.

Self-Tapping Bone Fixation Screws

1.5 mm x 8 mm

1 screw

PFB8

5 screws

PFB8-5

1.5 mm x 10 mm

1 screw

PFB10

5 screws

PFB10-5

1.5 mm x 12 mm

1 screw

PFB12

5 screws

PFB12-5

1.5 mm x 14 mm

1 screw

PFB14

5 screws

PFB14-5

 actual size

 actual size

 actual size

 actual size



Selection of Applicable References

Membrane References

- Wu IH, Bakhshalian N, Galaustian R, Naini RB, Min S, Freire M, Zadeh HH. Retrospective Analysis of the Outcome of Ridge Preservation with Anorganic Bovine Bone Mineral: Marginal Bone Level at Implants Placed Following Healing of Grafted Extraction Sockets. *Int J Periodontics Restorative Dent*. 2019 Jan/Feb;39(1):131-140.
- Phillips DJ, Swenson DT, Johnson TM. Buccal bone thickness adjacent to virtual dental implants following guided bone regeneration. *J Periodontol*. 2018 Dec 21.. [Epub ahead of print]
- Mendoza-Azpur G, Gallo P, Mayta-Tovalino F, Alva R, Valdivia E. A Case Series of Vertical Ridge Augmentation Using a Nonresorbable Membrane: A Multicenter Study. *Int J Periodontics Restorative Dent*. 2018 Nov/Dec;38(6):811-816.
- Sun DJ, Lim HC, Lee DW. Alveolar ridge preservation using an open membrane approach for sockets with bone deficiency: A randomized controlled clinical trial. *Clin Implant Dent Relat Res*. 2018 Nov 5.[Epub ahead of print]
- Changi KK, Greenstein G. Cytocone Procedure: Conservative Repair of a Buccal Plate Dehiscence in Preparation for Implant Placement. *Compend Contin Educ Dent*. 2018 May;39(5):294-299.
- Plonka AB, Urban IA, Wang HL. Decision Tree for Vertical Ridge Augmentation. *Int J Periodontics Restorative Dent*. 2018 Mar/Apr;38(2):269-275.
- Johnson M, Baron D. Tunnel Access for Guided Bone Regeneration in the Maxillary Anterior. *Clinical Advances in Periodontics*. Vol 8. No 1. March 2018.
- Bakhshalian N, Freire M, Min S, Wu I, Zadeh HH. Retrospective Analysis of the Outcome of Ridge Preservation with Anorganic Bovine Bone Minerals: Microcomputed Tomographic Assessment of Wound Healing in Grafted Extraction Sockets. *Int J Periodontics Restorative Dent*. 2018 Jan/Feb;38(1):103-111.
- Urban I, Traxler H, Romero-Bustillos M, Farkasdi S, Bartee B, Baksa G, Avila-Ortiz G. Effectiveness of Two Different Lingual Flap Advancing Techniques for Vertical Bone Augmentation in the Posterior Mandible: A Comparative, Split-Mouth Cadaver Study. *Int J Periodontics Restorative Dent*. 2018 Jan/Feb;38(1):35-40.
- Johnson TM, Berridge JP, Baron D. Protocol for Maintaining Alveolar Ridge Volume in Molar Immediate Implant Sites. *Clinical Advances in Periodontics*. November 2017, Vol. 7, No. 4, Pages 207-214
- Pistilli R, Checchi V, Sammartino G, Simion M, Felice P. Safe New Approach to the Lingual Flap Management in Mandibular Augmentation Procedures: The Digitoclastic Technique. *Implant Dent*. 2017 Oct;26(5):790-795.
- Urban IA, Monje A, Lozada J, Wang HL. Principles for Vertical Ridge Augmentation in the Atrophic Posterior Mandible: A Technical Review. *Int J Periodontics Restorative Dent*. 2017 Sep/Oct;37(5):639-645.
- Cucchi A, Vignudelli E, Napolitano A, Marchetti C, Corinaldesi G. Evaluation of complication rates and vertical bone gain after guided bone regeneration with non-resorbable membranes versus titanium meshes and resorbable membranes. A randomized clinical trial. *Clin Implant Dent Relat Res*. 2017 Jul 26. (epub ahead of print).
- Johnson TM, Berridge JP, Baron D. Protocol for Maintaining Alveolar Ridge Volume in Molar Immediate Implant Sites. *Clin Adv Periodontics*. 2017 Jul: 1-22.
- Gultekin BA, Cansiz E, Borahan MO. Clinical and 3-Dimensional Radiographic Evaluation of Autogenous Iliac Block Bone Grafting and Guided Bone Regeneration in Patients With Atrophic Maxilla. *J Oral Maxillofac Surg*. 2017 Apr;75(4):709-722.
- Ghensi P, Stablum W, Bettio E, Soldini MC, Tripi TR, Soldini C. Management of the exposure of a dense PTFE (d-PTFE) membrane in guided bone regeneration (GBR): a case report. *Oral Implantol (Rome)*. 2017 Nov 30;10(3):335-342.
- Laurito D, Lollobrigida M, Gianni F, Bosco S, Lamazza L, De Biase A. Alveolar Ridge Preservation with nc-HA and d-PTFE Membrane: A Clinical, Histologic, and Histomorphometric Study. *Int J Periodontics Restorative Dent*. 2017 Mar/Apr;37(2):283-290.
- Walker CJ, Prihoda TJ, Mealey BL, Lasho DJ, Noujeim M, Huynh-Ba G. Evaluation of healing at molar extraction sites with and without ridge preservation: a randomized controlled clinical trial. *J Periodontol*. 2017 Mar;88(3):241-249.
- Laurito D, Cugnetto R, Lollobrigida M, Guerra F, Vestri A, Gianni F, Bosco S, Lamazza L, De Biase A. Socket Preservation with d-PTFE Membrane: Histologic Analysis of the Newly Formed Matrix at Membrane Removal. *Int J Periodontics Restorative Dent*. 2016 Nov/Dec;36(6):877-883.
- Ronda M, Stacchi C. A Novel Approach for the Coronal Advancement of the Buccal Flap. *Int J Periodontics Restorative Dent*. 2015 Nov-Dec;35(6):795-801.
- Urban IA, Monje A, Wang HL. Vertical Ridge Augmentation and Soft Tissue Reconstruction of the Anterior Atrophic Maxillae: A Case Series. *Int J Periodontics Restorative Dent*. 2015 Sep-Oct;35(5):613-23.
- Al-Hezaimi K, Iezzi G, Rudek I, Al-Daafas A, Al-Hamdan K, Al-Rasheed A, Javed F, Piattelli A, Wang HL. Histomorphometric Analysis of Bone Regeneration Using a Dual Layer of Membranes (dPTFE Placed Over Collagen) in Fresh Extraction Sites: A Canine Model. *J Oral Implantol*. 2015 Apr;41(2):188-95.
- Borg TD, Mealey BL. Histologic healing following tooth extraction with ridge preservation using mineralized versus combined mineralized-demineralized freeze-dried bone allograft: a randomized controlled clinical trial. *J Periodontol*. 2015 Mar;86(3):348-55.
- Cucchi A, Ghensi P. Vertical Guided Bone Regeneration using Titanium-reinforced d-PTFE Membrane and Prehydrated Corticocancellous Bone Graft. *Open Dent J*. 2014 Nov 14;8:194-200.

Selection of Applicable References

- Ronda M, Rebaudi A, Torelli L, Stacchi C. Expanded vs. dense polytetrafluoroethylene membranes in vertical ridge augmentation around dental implants: a prospective randomized controlled clinical trial. *Clin Oral Implants Res.* 2014 Jul;25(7):859-66.
- Barboza EP, Stutz B, Mandarino D, Rodrigues DM, Ferreira VF. Evaluation of a dense polytetrafluoroethylene membrane to increase keratinized tissue: a randomized controlled clinical trial. *Implant Dent.* 2014 Jun;23(3):289-94.
- Urban IA, Lozada JL, Jovanovic SA, Nagursky H, Nagy K. Vertical Ridge Augmentation with Titanium-Reinforced, Dense-PTFE Membranes and a Combination of Particulated Autogenous Bone and Anorganic Bovine Bone-Derived Mineral: A Prospective Case Series in 19 Patients. *Int J Oral Maxillofac Implants.* 2014 Jan-Feb;29(1):185-93.
- Carbonell JM, Martin IS, Santos A, Pujol A, Sanz-Moliner JD, Nart J. High-density polytetrafluoroethylene membranes in guided bone and tissue regeneration procedures: a literature review. *Int J Oral Maxillofac Surg.* 2014 Jan;43(1):75-84.
- Vittorini Orgeas G, Clementini M, De Risi V, de Sanctis M. Surgical techniques for alveolar socket preservation: a systematic review. *Int J Oral Maxillofac Implants.* 2013 Jul-Aug;28(4):1049-61.
- Al-Hezaimi K, Rudek I, Al-Hamdan KS, Javed F, Nooh N, Wang HL. Efficacy of using a dual layer of membrane (dPTFE placed over collagen) for ridge preservation in fresh extraction sites: a micro-computed tomographic study in dogs. *Clin Oral Implants Res.* 2013 Oct;24(10):1152-7.
- Waasdorp, J, Feldman, S. Bone regeneration around immediate implants utilizing a dense polytetrafluoroethylene membrane without primary closure: A report of 3 cases. *J Oral Implantol.* 2013;39:355-361.
- Annibaldi S, Bignozzi I, Sammartino G, La Monaca G, Cristalli MP. Horizontal and Vertical Ridge Augmentation in Localized Alveolar Deficient Sites: A Retrospective Case Series. *Implant Dent.* 2012 Jun;21(3):175-185.
- Levin B. Immediate temporization of immediate implants in the esthetic zone: Evaluating survival and bone maintenance. *Compendium* 2011;32:52-62.
- Barboza EP, Stutz B, Ferreira VF, Carvalho W. Guided bone regeneration using nonexpanded polytetrafluoroethylene membranes in preparation for dental implant placements – A report of 420 cases. *Implant Dent.* 2010;19:2-7.
- Zafiropoulos GG, Deli G, Bartee BK, Hoffman O. Single-tooth implant placement and loading in fresh and regenerated extraction sockets. Five-year results: A case series using two different implant designs. *J Periodontol.* 2010;81:604-615.
- Zafiropoulos GG, Hoffmann O, Kasaj A, Wiltershausen B, Deli G, Tatakis DN. Mandibular molar root resection versus implant therapy: A retrospective nonrandomized study. *J Oral Implantol.* 2009;35:52-62.
- Fotek PD, Neiva RF, Wang HL. Comparison of dermal matrix and polytetrafluoroethylene membrane for socket bone augmentation: A clinical and histologic study. *J Periodontol.* 2009;80:776-785.
- Hoffman O, Bartee BK, Beaumont C, Kasaj A, Deli G, Zafiropoulos GG. Alveolar bone preservation in extraction sockets using non-resorbable dPTFE membranes: A retrospective non-randomized study. *J Periodontol.* 2008;79:1355-1369.
- Barber HD, Lignelli J, Smith BM, Bartee BK. Using a dense PTFE membrane without primary closure to achieve bone and tissue regeneration. *J Oral Maxillofac Surg.* 2007;65:748-752.
- Walters SP, Greenwell H, Hill M, Drisko C, Pickman K, Scheetz JP. Comparison of porous and non-porous teflon membranes plus a xenograft in the treatment of vertical osseous defects: A clinical reentry study. *J Periodontol.* 2003;74:1161-1168.
- Bartee BK. Extraction site reconstruction for alveolar ridge preservation. Part 1: Rationale and material selection. *J Oral Implantol.* 2001;27:187-193.
- Bartee BK. Extraction site reconstruction for alveolar ridge preservation. Part 2: Membrane-assisted surgical technique. *J Oral Implantol.* 2001;27:194-197.
- Lamb JW III, Greenwell H, Drisko C, Henderson RD, Scheetz JP, Rebitski G. A comparison of porous and non-porous teflon membranes plus demineralized freeze-dried bone allograft in the treatment of class II buccal/lingual furcation defects: A clinical reentry study. *J Periodontol.* 2001;72:1580-1587.
- Bartee BK. Evaluation of a new polytetrafluoroethylene guided tissue regeneration membrane in healing extraction sites. *Compend Contin Educ Dent* 1998;19:1256-1264.
- Bartee BK, Carr JA. Evaluation of a high-density polytetrafluoroethylene (n-PTFE) membrane as a barrier material to facilitate guided bone regeneration in the rat mandible. *J Oral Implantol.* 1995;21:88-95.
- Bartee BK. The use of high-density polytetrafluoroethylene membrane to treat osseous defects: Clinical reports. *Implant Dent.* 1995;4:21-26.

Suture References

- Silverstein LH, Kurtzman GM, Shatz PC. Suturing for optimal soft-tissue management. *J Oral Implantol.* 2009;35:82-90.
- Silverstein LH. Suturing principles: Preserving needle edges during dental suturing. *PPAD.* 2005;17:562-564.



INTERNATIONAL
EMAIL

4620 71st Street | Building 78-79
Lubbock, TX 79424
www.osteogenics.com
+1 806.796.1923
sales@osteogenics.com

Revised 3.19 EU | OBI CAT-002